
User Manual

CGM ESA-AD1101

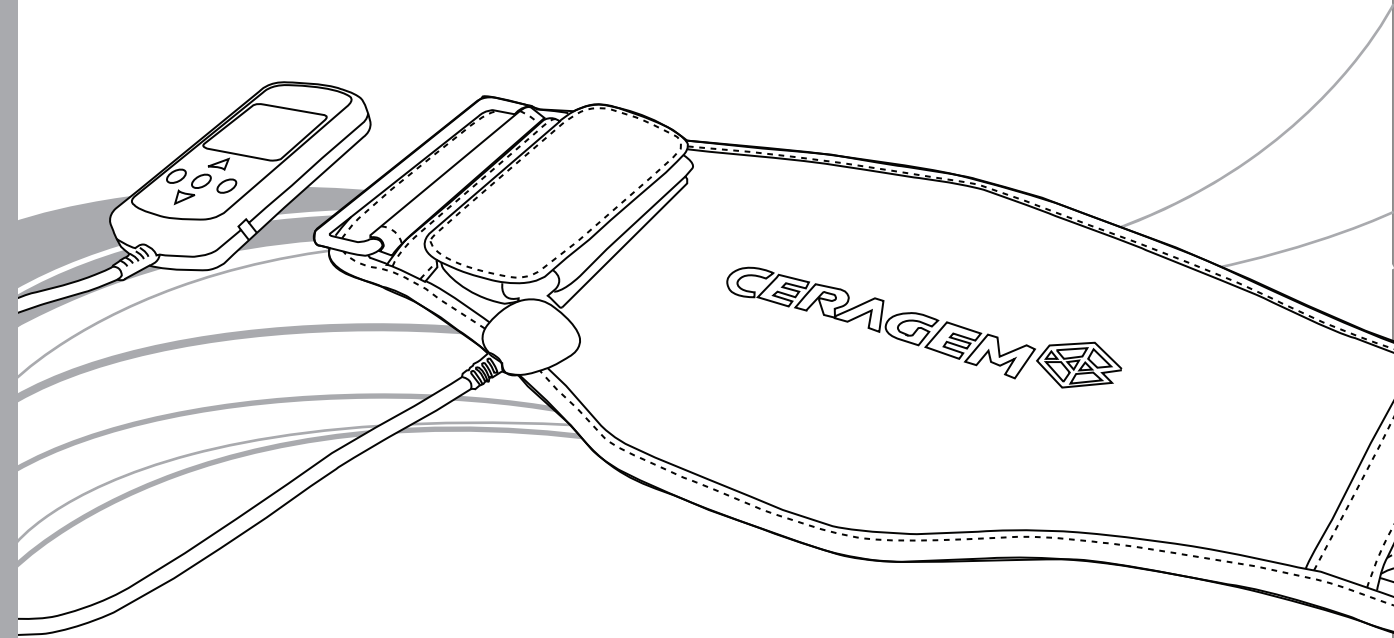Seller

CERAGEM CO., LTD.

177-14, Osaekdang-ri, Seonggeo-eup, Seobuk-gu,
Cheonan-si, Chungcheongnam-do, 331-831, KOREA
Tel : (82)41-523-2550 Fax : (82)41-523-2560
Homepage : <http://www.ceragem.net>

manufacturer

KUKDONG ELECTRONICS CO., LTD.

#411, 5Ra, Shihwa Industrial Complex, 673-17,
Seonggok-dong, Danwon-gu, Ansan-si, Gyeonggi-do, KOREA
Tel : (82)31-497-3356,3354 Fax : (82)31-497-3369

Please read this manual thoroughly for your safety.
Keep the manual in a safe place after reading it for future reference.


What is Slim Belt?


Slim Belt is a compact electronic physical massage machine for home use. Most today's people suffer from adult and chronic diseases due to irregular diet, various pollutions and stress. Though many of them try an exercise, dietary treatment, healthcare program or physiotherapy, it is hard to maintain it because of cost or space limitations. Slim Belt is a compact, multi-function electronic physical massage machine that you can use conveniently in your own home.

- Easy and convenient to use
- No danger of electricity and side effects
- No need to replace batteries thanks to the rechargeable battery.
- User-friendly directions on easy-to-read LCD screen
- Various massage effects with 3-channel massage functions.
- Best A/S service

Contents

1. Product Specifications and Components	2
2. Cautions for Using Slim Belt	3
3. How to Use Slim Belt's Controller	4
4. How to Use Controller Buttons	6
5. How to Use Rechargeable Battery	7
6. Precautions for Safe Use	8
7. Home Maintenance	9
8. Product Warranty	10

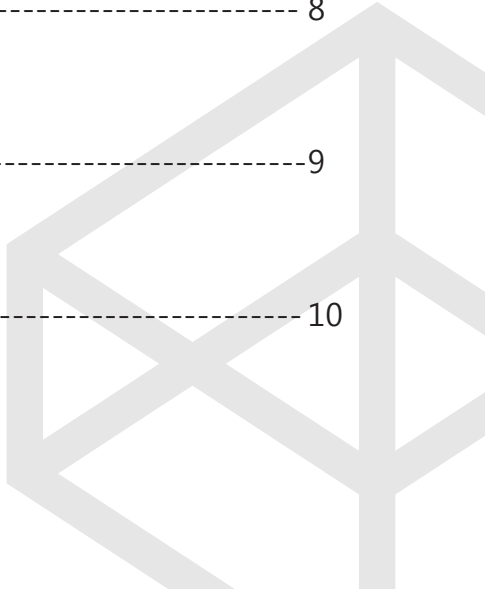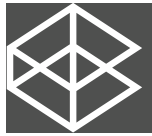1. Product Specifications and Components

Product Name	Slim Belt
Model	CGM ESA-AD1101
Display	Liquid crystal display (39 x 40 mm)
Frequency	1~200Hz
Application	Constant voltage type
Rated Voltage	DC3.7V (Li-Ion rechargeable battery)
Classification	Embedded power type
Operation Time	5~60 minutes
Operation Mode	1~6 levels
Polarity Reversal	Automatic polarity reversal LCD
Size	Controller (56 x 118 x 19 mm) Slim Belt (490 x 200 x 10mm)
Weight	1.1 kg (including box)
Components	Slim Belt (including controller) 1ea Pouch 1ea Recharging adaptor 1ea User's manual/warranty 1 copy

2. Cautions for Using Slim Belt

When to Use

We recommend using Slim Belt between meals or (more than 1 hour) after a meal but not on an empty stomach. Slim Belt is most effective after a bath.

Cautions for Using

Application of electric current to the same location for extended periods may cause fatigue. Please note that use of electric current stronger than recommended current settings may reduce its effectiveness.

Period to Use

For optical effect, use Slim Belt for 20~30 minutes 1~2 times daily. Reduce the time of use by 1/2 for beginners or 1/3 if you are suffering from low blood pressure, anemia or in a weak condition.

Overuse, especially during initial application, may cause fatigue or temporary elevation of body temperature.

For Minors (Children)


Minors (children) may use it but at lower settings and for shorter times.

- Under 3: Not recommended.
- 3~6 years old: Adjust intensity and time to 1/2 of recommendations for adults.
- 7~12 years old: Adjust intensity and time to 2/3 of recommendations for adults.
- Over 13: Follow recommended intensity and time for adults.


3. How to Use Slim Belt's Controller

Controller

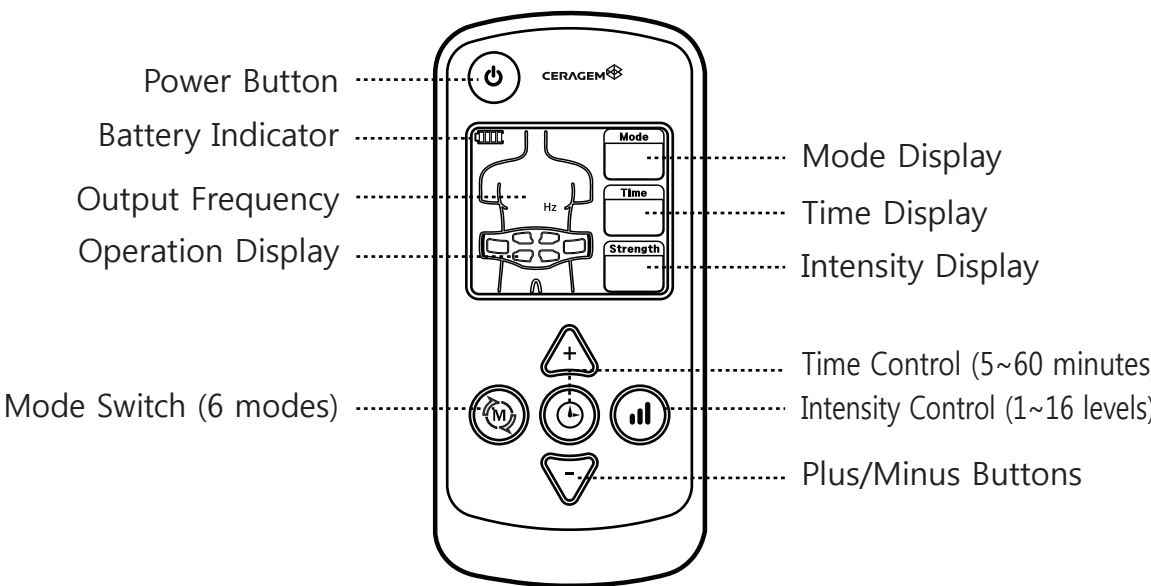Ergonomically Designed Program and Pad

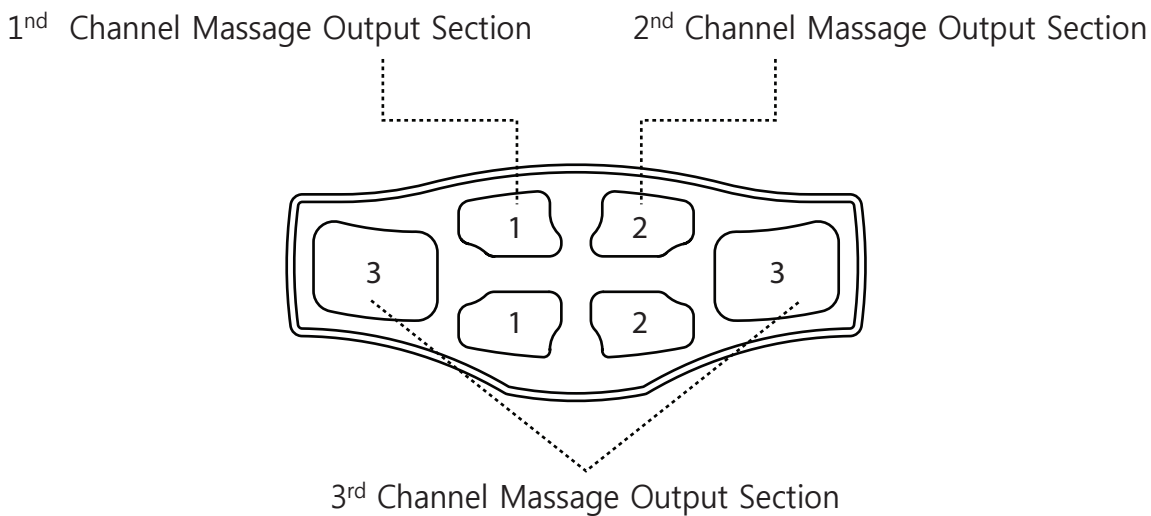3. How to Use Slim Belt's Controller

Power Button

- When the power is off, press and hold the power button (for more than 1 second); it will beep once and the power will be turned on.
- When the power is on, press and hold the power button (for more than 1 second); it will beep twice and the power will be turned off.

Battery Indicator

- Indicates current battery status.

Time Display

- Indicates remaining operation time.

Output Frequency

- Indicates current frequency.

Mode Display

- Indicates current mode.

Operation Display

- Indicates location of operating pad (channel).

Intensity Display

- Indicates current output intensity.

4. How to Use Controller Buttons

1. Time Control (5~60 minutes)
1. Press the time button () once; the product will beep and the current time will flicker for 5 seconds.

2. You can change time by pressing plus/minus buttons ( ).

3. The product will be turned off automatically after the selected time.
2. Intensity Control (1~16 levels)
1. Press the intensity button () once; the product will beep and the current intensity level will flicker for 5 seconds.

2. You can change output intensity to 1~16 levels by pressing plus/minus buttons ( ).
3. Mode Switch (6 modes)
1. Press the mode button () once; the product will beep and the current mode number will flicker for 5 seconds.

2. You can change MODE to 1~6 by pressing plus/minus buttons ( ).
(Pressing MODE button repeatedly also change the mode.)

Mode 1	Cycle mode 2~6	Mode 2	Fast beat
Mode 3	Slow beat	Mode 4	Beat
Mode 5	Fast beat 2	Mode 6	Massage

5. How to Use Rechargeable Battery

1. Specifications of rechargeable battery
- Battery : Li-Ion


Operating Voltage : 3.7 V

Recharging Capacity: 1,200 mA

Recharging Voltage : 4.2 V
2. Recharging adaptor
- Input Voltage : AC100-240 V

Output Voltage : 4.2 V

Output Current : 1,000 mA
3. How to recharge
- If the power of Slim Belt is insufficient or turned off, recharge it as the following diagram.

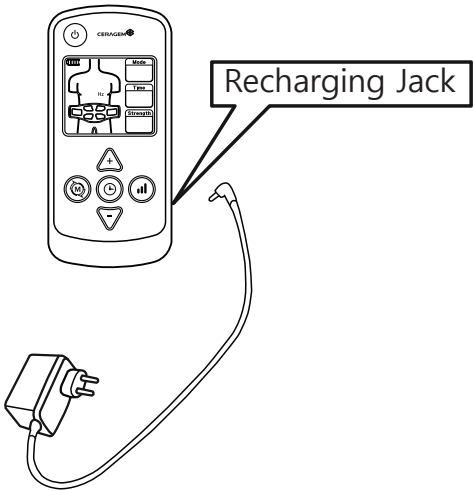- Connect the recharging adaptor to the side of adjustor. The LED lamp of the recharging adaptor will flicker in red.
- When fully charged, the LED lamp will flicker in green.
- [Recharging Time: about 3 hours]

6. Precautions for Safe Use

Caution

The cautions contained in this manual indicate a potentially hazardous situation that may cause minor injury or property damage.

Caution


1. Do not use this product with the following medical electronic equipment.
 - Medical equipment inserted into body such as pacemaker
 - Life support equipment such as heart-lung machine
 - Mounted type medical equipment such as cardiograph
 - Other similar medical electronic equipment.
2. Consult your doctor before using Slim Belt if you:
 - have a heart disease / acute disease / malignant tumor / high fever
 - have a perception disorder / are pregnant / have a skin disease where Slim Belt will be applied.
 - have been recommended by a doctor not to use this product.


1. When not in use, replace in the original container and store in a safe place.

2. Do not attempt to dismantle, repair or modify this product.

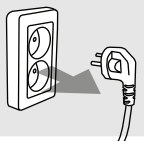3. Always turn off the power after use.

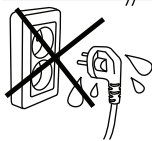4. Do not touch the power cord with wet hands.

5. Do not wash the product with benzene, thinner or acetone.

6. Prevent children from sucking or playing with this product when it is turned on.

7. Home Maintenance

1. LCD is not turned on.

Battery is discharged or not connected properly.

- Recharge or replace the battery.

2. Failure to feel stimulus.

Battery is discharged or adhesive pad is dead.

- Apply the resolution specified in No. 1.

3. Stimulus is weak or the product stings.

Adhesion surface of pad is dirty or contaminated severely.

- Clean the adhesive pad with a wet towel before use. Adhesive pad is not attached to the body completely.
- Spread massage gel or water evenly on the surface of adhesive pad.
- Please call A/S center

4. Stimulus is weak or unfelt even at a high intensity.

- Clean the adhesive pad with a wet towel before use.
- Apply the resolution specified in No. 3.

5. Power becomes weaker after 1 month.

Battery becomes weak.

- Apply the resolution specified in No. 1.

Your body becomes accustomed to the stimulus.

- Increase the intensity.

8. Product Warranty

Warranty

Product Name: Slim Belt
Model: CGM ESA-AD1101

Year/Month/Day of Purchase			
Place of Purchase			
Warranty Period		1 year from the date of purchase	
Customer Information	Name	Contact	
	Address		

This product is manufactured and tested by thorough product quality management standards.
If you experience any spontaneous failure, contact the dealer or our A/S team with this warranty.

Caution

- Bring this warranty because without it, the replacement or repair of the product will be impossible.
- Complete the warranty form on the date of purchase.
- For other questions, please contact the dealer or our A/S center.
- To improve quality, design and product specifications may be changed without notice.

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